



TR 9650

Mobile bath lift system

Hygiene Equipment
Designed for People

TR 9650 Designed for maximum flexibility

The TR9650 is a mobile bath lift system designed to easily convert into different configurations, to accommodate a wide variety of needs for patients and residents in both the acute and long term care environments.

Whether in standing, seated or recumbent position, the TR 9650 provides a dignified complete solution for patients needing aid getting in or out of the bathtub.

Easy access and comfortable padded seating surfaces ensure a pleasant bathing experience.

Ergonomically designed hand grips and safetybelts are incorporated on the product, and provides security and assurance for both the patient and the caregiver.

TR9650 can be converted between a chair and a stretcher without the need of any tools. This makes the product a good choice for multiple mobility levels.



Registered at Medical Products Agency, Sweden as a MDR Class 1 product according to MDR 2017/745.

TR reserves the right to make design changes and/or improvements without notice.





TR 9650

Mobile bath lift system

Hygiene Equipment
Designed for People

- Foot operated and height adjustable enables the caregiver to find the most ergonomic working height and also gives a hands free opportunity to tend the patient.
- Seating and stretcher surfaces are made of a soft and durable water absorbent material. This makes it easy to clean, maintain and is comfortable to sit and lay on.
- Removable side parts provide a simple configuration from chair to stretcher, with 3 angles, 27, 40 or 54 degrees.
- The safety arm support in both the chair and stretcher configurations can be folded away for easy access and moved back while the patient is in position.
- The leg configuration is designed for easy access under beds and baths.
- An integral autostop ensures that the seat or stretcher cannot drop if the lifting belt breaks.
- TR9650 can be used with almost all supine bath systems. An optional spacer can be fitted in case needed.

Specifications


Operating height	200–965 mm
Working load limit	150 kg
Castors	4 twin castors, 100 mm as standard. 2 locking and 2 straight steering
Control	Foot operated control for height.
Battery	Rechargeable batteries 24 V 7Ah

Models:

- TR 9650 Chair
- TR 9650 Stretcher
- TR 9650 Combi model
- TR 9650 Side parts for combi model

Options/Accessories:

- Removable bed pan holder and bucket
- Extension adapter to increase distance from lift column to chair/stretcher

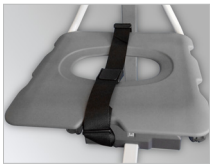


Opening for bucket



Easy to clean



Charger/emergency stop

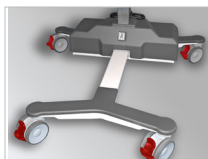


Chassi




Side parts

Dimensions (mm)

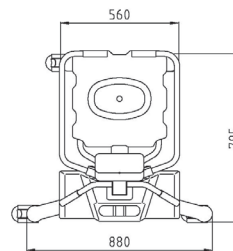

Chairlift

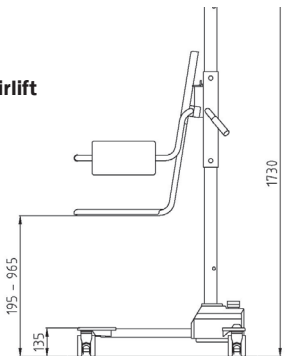

Stretcher/Combilift

